

E-track

Overhead tracking: Go to the wall

Primary application

Home care
Nursing homes
Hospital

Ergolet's E-Track is a proud example of innovative Danish design. It represents, quite simply, a revolution in overhead tracking techniques.

E-Track's unique difference is that it is a wall mounted track solution. In reality E-Track is an attractive flat wall panel offering a discreet and unobtrusive alternative to the unsightly wall and ceiling brackets that are conventionally used. Moreover, E-track is quick and simple to install, making it an extremely cost effective choice and the preferred professionals' choice for home and institutional settings.

E-Track is also a supremely flexible system. For example, a room can be completed with E-Track and a traverse system, but if the traverse system is not being used, or is required in another room, then the traverse can be quickly and easily dismantled from the E-track profiles and conveniently moved.

An additional feature is E-Track's profile on the underside, a practical feature for hanging light pictures or other wall decorations, removing the need to make holes in the wall.

An aluminum cover plate locks onto the lower fixing channel concealing unsightly electrical cables and screws. And if a room has E-Track installed but is not being used for overhead lifting, then a second cover plate can be clicked into the upper channel transforming E-Track it into a flat, fully closed wall panel.

Key standard functions


Can be mounted on almost any type of wall construction

Accessories

A range of accessories are available for E-track such as brackets for hanging pictures etc.

A range of brackets for different installation purposes.




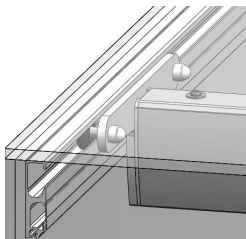


Unique telescopic cart to compensate for varying wall widths.




Quick and easy to install. Saves time and reduces costs.

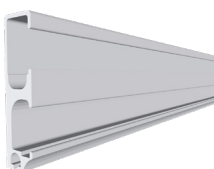


Discreet thin profile.

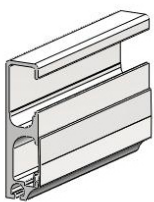

The profile has an extra track underneath, which can be used to hang pictures, book shelves or other wall decorations.

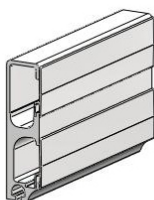

DOM
OVER

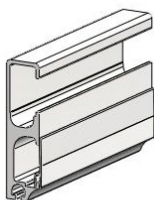

EEC




PRODUCT HIGHLIGHTS

- Class leading design and functionality.
- Can be mounted on most wall constructions.
- Quick and easy to install. Cost saving.
- Telescopic cart to compensate for varying wall widths.
- Profile can be used to support light room dividers or curtains.
- Clean and simple state-of-the-art-design. Flat panel look.

Version	Item-number
E-track 4,3m	016-01028
E-track 5,5m	016-01029
E-track 7m	016-01027
E-track 7,5m	016-01026

TECHNICAL SPECIFICATION

Load on profile	200kg and 275kg
Measurements of profile	25,5 mm x 97 mm
Weight	2,0 kg/m
Mounting examples	Plasterboard, wood, concrete, brick
Colour	White powder coating

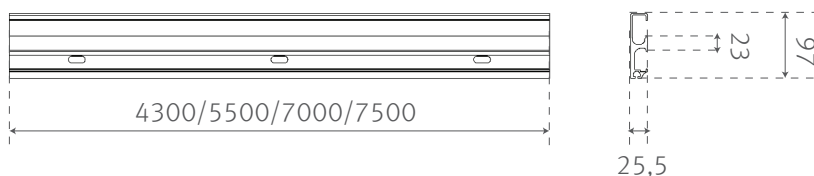

Consult the Technical manual for further specifications. All figures and measurements are for guidance only and may differ in actual situations. Ergolet reserves the right to changes without notice.

Product photos may include accessories. Consult your dealer.

WinnCare Nordic ApS
 Taarnborgvej 12C
 4220 Korsoer
 Denmark
 tel. +45 70 27 37 20
 e-mail: info@winnCare.dk

WINNCARE
 N o r d i c

www.winnCare.dk