


ALOVA™ DUO



Braden scale score

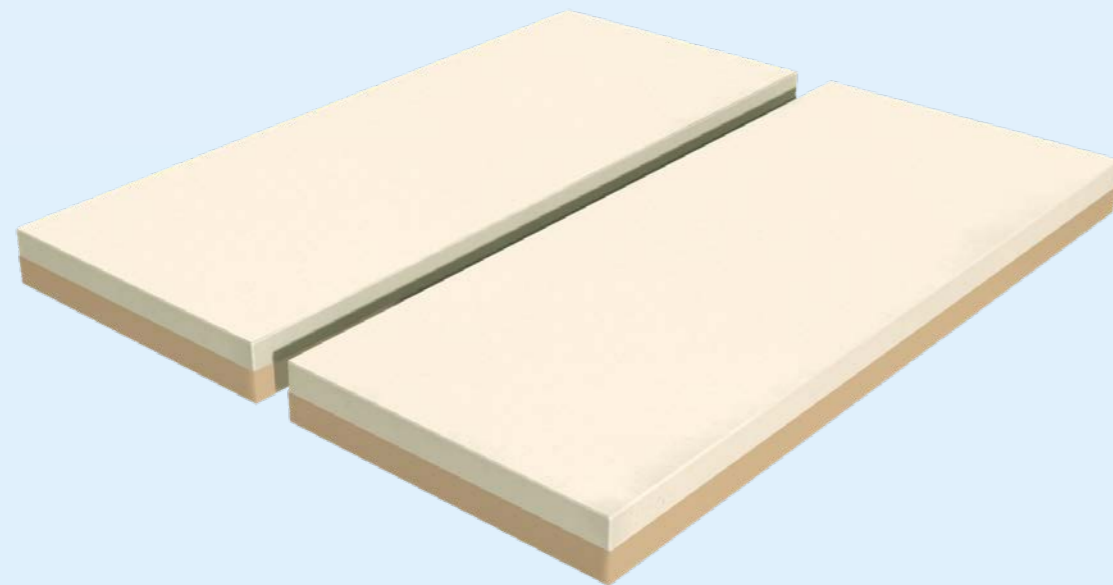

KEY FEATURES

- ▶ Single-piece mattress made up of:
 - 7 cm molded visco elastic top layer: density: $\geq 75 \text{ Kg/m}^3$; hardness: 2.2 kPa
 - 7 cm High Resilient (HR) foam base: density: 34 Kg/m^3 ; hardness: 1.8 kPa
- ▶ Patient weight: 30 to 135 Kg
- ▶ Dimensions: VMA70: 196 x 70 x 14 cm (77.1 x 25.5 x 5.5 inches)
VMA80: 196 x 80 x 14 cm (77.1 x 31.4 x 5.5 inches)
- ▶ Mattress weight: 16.8 Kg (VMA70) and 18.2 Kg (VMA80)
- ▶ Cover available: Promust PU HD

KEY BENEFITS

- ▶ Visco elastic foam provides excellent pressure re-distribution properties: homogeneous immersion of the body into the foam, thus increasing the bearing surface and reducing pressure applied in areas at risk of pressure ulcers .
- ▶ Suitable for individual sleeping surfaces of 140 or 160 cm width medical beds for couples hospitalized at home.
- ▶ Visco elastic foam provides patient comfort.

INDICATIONS

HELP IN PREVENTING PRESSURE ULCERS FOR PATIENT AT MEDIUM TO HIGH RISK OF PRESSURE ULCER DEVELOPMENT, BEDRIDDEN MORE THAN 15 HOURS A DAY.

HELPS INTREATING EXISTING PRESSURE ULCERS FROM STADE 1 TO 2 COMBINED WITH TECHNICAL AID DEVICES.

PATIENT WEIGHT: 30 TO 135 KG/ 4.7 TO 21.2 ST

- ▶ Single-piece mattress made up of visco elastic foam top layer on High Resilient (HR) polyurethane foam
- ▶ ALOVA DUO: therapy and comfort in bed for couples
- ▶ Removable and washable integral cover

DESCRIPTION	ITEM REFERENCES ACCORDING TO PACKAGING CARDBOX	FIRE STANDARDS
ALOVA™ DUO 70 cm for 140 cm wide medical bed with double bed frame, integral cover in Promust PU HD	VMA70	EN 597-1 and 2 GPEM D1 90 classe D
ALOVA™ DUO 80 cm for 160 cm wide medical bed with double bed frame, integral cover in Promust PU HD	VMA80	

3 years warranty

Cleaning instructions: Promust PU HD cover

FP000037/01.2018

Photos not contractually binding

